

INTUY[®] KNEE

PRODUCT HIGHLIGHTS:

- POWERFUL, LIGHTWEIGHT
- IP45 WEATHERPROOF
- LOW NOISE
- 44HR BATTERY STANDBY TIME
- 2HR FULL CHARGE
- REGENERATIVE CHARGING
- 100NM TORQUE
- EFFORTLESS MODE TRANSITION
- SMOOTH, NATURAL AMBULATION

OVERVIEW:

The Intuy[®] Knee by Rebocon is designed to promote natural gait biomechanics and reduce strain on sound joints by employing advanced processors, sensors and algorithms.

The lightest, quietest, and most powerful prosthetic knee available, the Intuy Knee uses intuitive symmetrical sitting, a standing positional lock, and powered extension, to deliver consistent natural gait and stability to a wide user range with its 100Nm torque, 275lb weight limit, and 255mm build height. The knee also features seamless transitions into stair ascent/descent, sit-to-stand, and biking assistance functions.

44-hour autonomy, a full charge from empty in just 2 hours, combined with regenerative charging during resistance encounters (sitting, descending stairs or ramps) provides the patient with reliable, multi-day support.

Product is only available direct through WillowWood.
Contact intuyorders@willowwood.com for more information.

INTUY KNEE: SPECIFICATIONS

WEIGHT LIMIT: 124.74kg (275lbs)

BUILD HEIGHT: 255mm (10.04")

PRODUCT WEIGHT: 2.4kg (5.29lbs)


*w/ battery

WARRANTY: 3 years
from date of patient fitting

KNEE FLEXION: 120°



PDAC APPROVED L-CODES:
L5828, L5845, L5848, L5856, and L5859

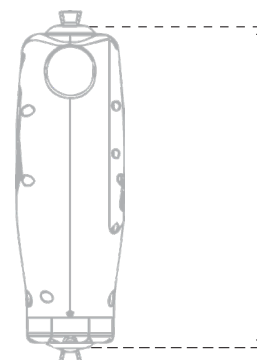

INTUY KNEE
BUILD HEIGHT
255 MM (10.04")

WillowWood

PN-2717-A ©WillowWood Global LLC, 2023