

ALPHA CONTROL™ LINER SYSTEM

PRODUCT HIGHLIGHTS:

- IMPROVES USER CONTROL AND ALLOWS 11% GREATER ELBOW RANGE OF MOTION
- ELIMINATES THE COMPLEX PROCESS OF ELECTRODE PLACEMENT
- IMPROVES SIGNAL CONSISTENCY BY MINIMIZING ELECTRODE LIFT OFF
- IMPROVED SOCKET COMFORT AND EASIER DONNING & DOFFING FOR PATIENTS
- EASES FABRICATION PROCESS
- FOR USE WITH COAPT GEN2 COMPLETE CONTROL SYSTEM



OVERVIEW

The Alpha Control Liner System is the next generation interface solution for use with pattern recognition. Embedded electrodes and leads in a gel-covered locking liner simplify the process of connecting to the Coapt Gen2 Complete Control System. The liner is available in three sizes (Extra Small, Small and Medium), with nine potential configurations for lead length – ranging from 6 cm to 22 cm (in 2 cm increments) – enabling simplified myo-site selection and accommodate the vast majority of residual limb sizes.

By enabling consistent signal capture (through elimination of electrode lift-off), the system delivers a hassle-free user experience – enabling users to expand their range of activities while ensuring superior intuitive control and comfort.


The system is easier to don and doff in comparison with standard skin fit interface methods – leading to higher user satisfaction with their prosthetic device.

The Alpha Control Liner System was developed in collaboration with Coapt and Shirley Ryan AbilityLab and is designed for the external prosthetic fittings of upper limbs and functions exclusively with the Coapt Gen2 Complete Control System.

- ✔ TRANSRADIAL
- ✔ TRANSHUMERAL

LCODE SUGGESTIONS: L7499

Liner Sizing:

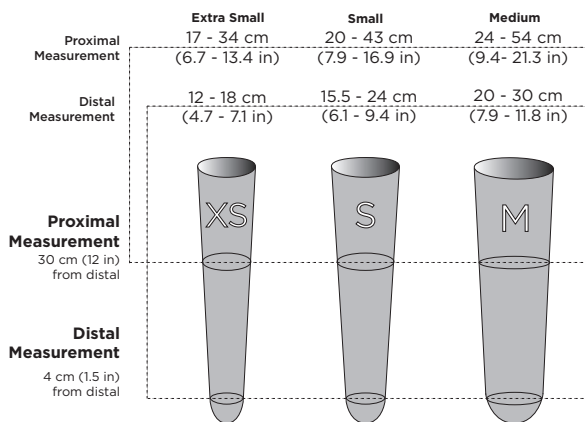

ALPHA CONTROL SYSTEM PART NUMBERS

Leads	Extra Small	Small	Medium
6 cm	ACLS-049-X06	ACLS-049-S06	ACLS-049-M06
8 cm	ACLS-049-X08	ACLS-049-S08	ACLS-049-M08
10 cm	ACLS-049-X10	ACLS-049-S10	ACLS-049-M10
12 cm	ACLS-049-X12	ACLS-049-S12	ACLS-049-M12
14 cm	ACLS-049-X14	ACLS-049-S14	ACLS-049-M14
16 cm	ACLS-049-X16	ACLS-049-S16	ACLS-049-M16
18 cm	ACLS-049-X18	ACLS-049-S18	ACLS-049-M18
20 cm	ACLS-049-X20	ACLS-049-S20	ACLS-049-M20
22 cm	ACLS-049-X22	ACLS-049-S22	ACLS-049-M22

*Liner is only available in a 3 mm uniform profile

FOR MORE INFORMATION CALL 800.848.4930 OR 740.869.3377
OR VISIT WILLOWOOD.COM