

ALPHA TRUSEAL LINER

PRODUCT HIGHLIGHTS:

- MEDICAL GRADE, PLATINUM-CURED SILICONE
- FOUR INTEGRATED SEALING FINS
- SATIN OUTER-FINISH, NO LUBRICANTS NEEDED
- EMBEDDED FABRIC MESH
- FEATURES A 12-MONTH WARRANTY



Alpha Silicone VMP, set of 2
Pad measures 9" x 4.5"; 4.5 mm & 7.5 mm thicknesses
ALA-VM-S2


OVERVIEW

WillowWood's Alpha TruSeal Liner uses a proprietary, platinum-cured, medical grade silicone designed for comfort and performance. The liner has four integrated sealing fins designed to ensure confident suction suspension. With TruSeal's stepped fins and Volume Management Pads, a large range of fluctuation can be accommodated. The satin outer finish allows for easy donning, does not require messy lubricants and is easier to clean than a fabric covered liner. An inner topcoat is applied to reduce tack while controlling soft tissue compression and rotation. An embedded mesh fabric within the silicone resists abrasion and protects the liner from rips and tears. The Alpha TruSeal Liner is available in a symmetrical profile and two tapers (conical and cylindrical) to accommodate various limb shapes.

✓ TRANSFEMORAL

LCODE SUGGESTIONS: L5679 (cushion)


Liner Sizing: Conical

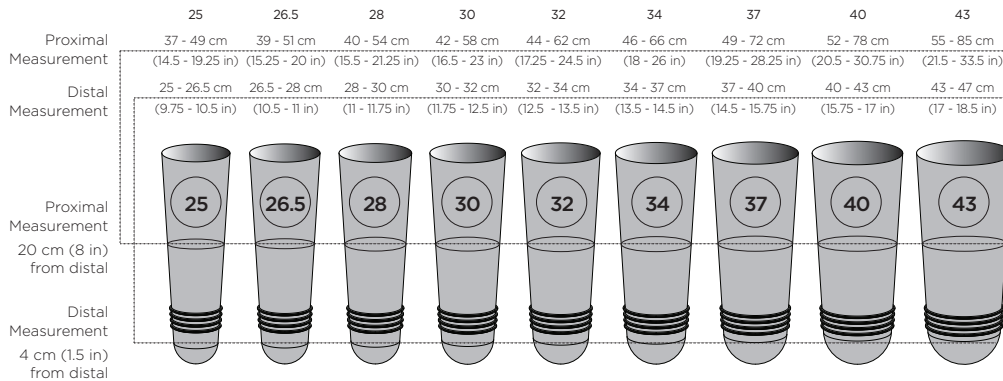

Size	Part #
25	S494-1125
26.5	S494-1126
28	S494-1128
30	S494-1130
32	S494-1132
34	S494-1134
37	S494-1137
40	S494-1140
43	S494-1143

Optional
ALA-VM-S2 Alpha Silicone Volume Management Pads, set of 2

Liner Sizing: Cylindrical

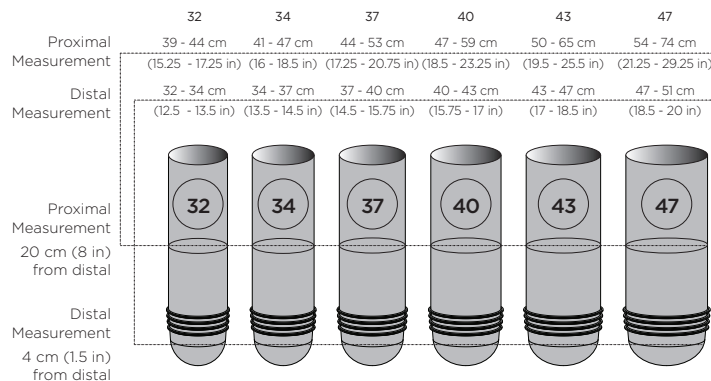

Size	Part #
32	S495-1132
34	S495-1134
37	S495-1137
40	S495-1140
43	S495-1143
47	S495-1147

Optional
ALA-VM-S2 Alpha Silicone Volume Management Pads, set of 2