

WILLOWWOOD®

INTUY® Knee



Product Highlights

- Powerful, lightweight
- IP45 weatherproof
- Low noise
- 44 hr battery life (standby time)
- 2 hr full charge
- Regenerative charging
- 100 Nm torque
- Effortless mode transition
- Smooth, natural ambulation

The INTUY Knee by Rebocon is designed to promote natural gait biomechanics and reduce strain on sound joints by employing advanced processors, sensors, and algorithms.

The lightest, quietest, and most powerful prosthetic knee available, the INTUY Knee uses intuitive symmetrical sitting, a standing positional lock and powered extension, to deliver consistent natural gait and stability to a wide user range with its 100 Nm torque, 275 lb weight limit, and 255 mm build height. The knee also features seamless transitions into stair ascent/descent, sit-to-stand, and biking assistance functions.

Up to 2 days of battery autonomy, a full charge from empty in just 2 hours, combined with regenerative charging during resistance encounters (sitting, descending stairs, or ramps) provides the patient with reliable, multi-day support.

Product is only available direct through WillowWood. Contact intuyorders@willowwood.com for more information.



PDAC Approved L-Codes:
L5828, L5845, L5848, L5856, and L5859

INTUY App



INTUY Knee: Specifications


Weight Limit: 124.74 kg (275 lbs)

Build Height: 255 mm (10.04")

Product Weight: 2.4 kg (5.29 lbs) *w/ battery

Warranty: 3 years from date of patient fitting

Knee Flexion: 120°

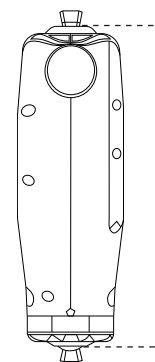

INTUY Knee
Build Height
255 mm (10.04")

FOR MORE INFORMATION CALL 800.848.4930 OR 740.869.3377
OR VISIT WILLOWWOOD.COM

<https://www.willowwood.com/education-resources/resources/patents/>

PN-2717-E

WILLOWWOOD.

©WillowWood Global LLC, 2024