

WILLOWWOOD®

META® Flow™

Product Highlights



- Polycentric ankle with carbon composite ankle spring – guaranteed freedom from fluid leaks
- 95% energy efficient – only design that eliminates dampening effect, provides power to propel user through mid-stance
- Delivers 14° Ankle RoM (9° plantarflexion + 5° dorsiflexion)
- Significantly reduces fall risk via 5° active dorsiflexion during swing phase

WillowWood's META® Flow™ is the first PDAC approved L5968 foot that eliminates the trade-off between range of motion and dynamics. The design combines the benefits of the META Unibody® platform with an innovative polycentric ankle featuring a carbon composite ankle spring. The unique design is the only L5968 foot that provides energy during midstance – allowing users to experience best-in-class rollover and freedom from the dampening effect experienced with hydraulic designs. The end result is a smooth natural gait on various terrains and at varying cadences while reducing the risk of falls through active swing phase dorsiflexion. By eliminating the need for hydraulic fluid, the design ensures freedom from leaks and reduces the need for frequent maintenance.

The META Flow presents the optimal solution for K3 users across all impact levels seeking ankle mobility and dynamic response tailored to a diverse range of activities.



**PDAC
APPROVED**



- ✓ Transfemoral
- ✓ Transtibial

PDAC Approved L-Code: L5981, L5968

META Flow Foot: Specifications

Heel Height: 10 mm (3/8")


Build Height: 112 mm (4.41")

Approximate Product Weight: 862 g

*Std. size 26 Cat. 5 with foot shell and sock

Trial Period: 60 days from date of patient fitting

Warranty Period: 36 months from date of patient fitting

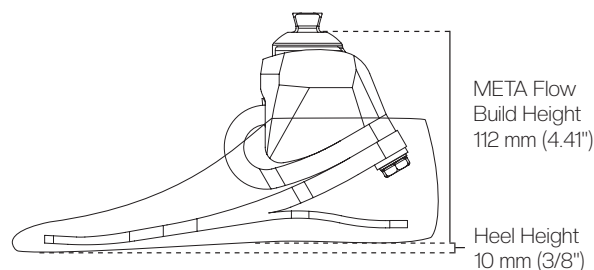

META Flow
Build Height
112 mm (4.41")

Heel Height
10 mm (3/8")

Sandal Toe Foot Shells



Medium Brown



Tan



Buff

META Flow Category Chart

Adjusted User Weight (patient body weight + loads carried)	Standard (Low to Medium) Impact Level		High Impact Level	
	Sizes 22-24 cm	Sizes 25-30 cm	Sizes 22-24 cm	Sizes 25-30 cm
99-115 lbs (45-51 kg)	1	1	1	1
116-130 lbs (52-58 kg)	1	1	2	2
131-150 lbs (59-67 kg)	2	2	3	3
151-170 lbs (68-76 kg)	3	3	4	4
171-194 lbs (77-87 kg)	4	4	5	5
195-221 lbs (88-99 kg)	5	5	-	6
222-256 lbs (100-115 kg)	-	6	-	7
257-287 lbs (116-130 kg)	-	7	-	8
288-324 lbs (131-147 kg)	-	8	-	9
325-365 lbs (148-165 kg)	-	9	-	-

*Note: The Standard impact level is recommended for most patients for day-to-day use; high impact level should be chosen for active outdoor users engaged in activities that involve heavy ground impact forces (e.g. sports, construction work)

The META Lineup

All META feet feature the groundbreaking META-Unibody® design that combines dynamic response with lightweight, hardware-free design suitable for K3-K4 users. The design delivers a smooth roll-over experience and delivers support in mid-terminal stance enhancing stability without compromising function.



META Core

Smooth rollover and late-stance support for K3-K4 users low-build height and weight.

Available in Carbon Fiber or Fiberglass.*

PDAC approved: L5981



META Arc

Biomimetic ankle offers stability and adaptation to terrain. 20° of inversion & eversion.

PDAC approved: L5981 + L5986



META Shock

Adaptive shock absorption via integrated shock module 13 mm of vertical displacement.

PDAC approved: L5987



META Shock X

Adaptive shock absorption and impact protection via 40° of torsional rotation + 13 mm vertical displacement.

Available in Carbon Fiber or Fiberglass.*

PDAC approved: L5987 + L5984



META Flow

Active dorsiflexion foot enables natural and safe gait over uneven terrain. 9° of plantarflexion & 5° of dorsiflexion.


PDAC approved: L5981 + L5968

How To Order

META Flow Foot with Spectra Sock

Ordering example:

META Flow Foot, 25 cm, left foot, toe resistance 5

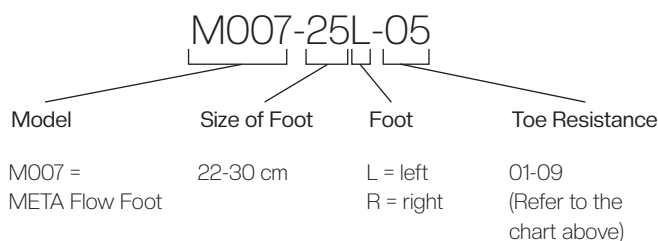

META Flow Foot Shell

Ordering example:

META Flow Foot Shell, 25 cm, left foot, tan

